TEST HIV

Test destiné à l'auto-diagnostic

Exacto

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TEST HIV
For self-testing

EXACTO® TEST HIV is a screening self-test for the detection of an infection by the HIV virus (HIV-1 or HIV-2), the virus responsible for AIDS (Acquired Immunodeficiency Syndrome). The test contains antigens capable of detecting antibodies directed against the virus responsible for AIDS (HIV-1 or HIV-2). It must be used at least three months after the last risky exposure to give a reliable result.

CONTENT OF THE KIT

Ref 855164 contains one bag for one self-test.
Ref 855169 contains two bags, one self-test can be run with each bag.

ADVICE FOR SAMPLING

First, wash your hands with hot water and dry them properly. Warm water helps improve blood circulation at the fingertips. You can let your hand loose along your body for around 1 minute in order for blood to flow in that area in larger quantities. It is recommended to puncture the left hand for a right-handed person, and the right hand for a left-handed person. Please puncture preferably a side of the fingertip: it is where the blood flow is higher, and the sensitivity minimal. Avoid areas where the skin seems thicker, so more difficult to puncture. It is possible that a blood drop does not appear immediately. In that case, please wait a little bit and massage lightly the finger from its base to the punctured area. If the red square is not completely filled with blood (step 12), you can form a new drop of blood, collecting it with the same inverted cup, and then deposit this new drop in the square «BLOOD».

Should you have any questions on how to run and interpret the test or want further information on the virus responsible for AIDS, its modes of transmission of the virus or screening centres :

WWW.SIDA-INFO-SERVICE.ORG
SIDA INFO SERVICE : +33 (0) 800 840 840 (24/7)

BLOOD COLLECTION

Visit www.exacto.fr to watch an instructional video or give your feedback about EXACTO® Test HIV

- 1** Please read the instructions thoroughly. Have a watch or a timer (only equipment not provided).
- 2** Wash your hands with soap and hot water and then dry them before proceeding to the next step.

Take all the elements out of the box, place them on a flat and clean surface and identify each element of the kit (see content on the left, 7 elements to count).


Open the pouch **A** and pull the cassette **B** out. Use the test within 1 hour

Pull the alcohol wipe **E** out of its pouch; and disinfect the tip of the finger chosen for blood collection.

Take the sterile pad **F** out of its sachet. Wipe away traces of alcohol.

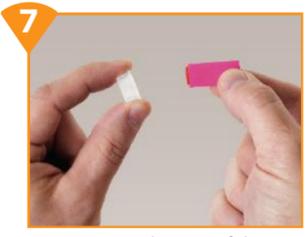
Remove the cap of the lancet **G**.

Apply the lancet **G** on the fingertip. Trigger it by pressing firmly on the other end.

Squeeze the tip of the finger to form a **large drop of blood** for the sampling.

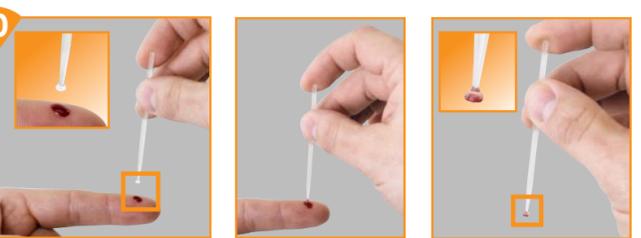Place the flared part of the inverted cup **D** on the drop of blood to fill entirely the cup.

TEST PROCEDURE

Immediately deposit the drop of blood into the square well BLOOD of the cassette **B** (The blood drop must be in direct contact with the bottom of the well).

Check that the bottom of the square well BLOOD is completely red before moving to the next step (this confirms that the drop of blood has been properly added).


Unscrew the cap of the dropper bottle **C**. Put two drops in the round well DILUENT of the cassette **B**.


Start the timer or write down the time of deposit. Wait 10 minutes to read the result. Do not read beyond 20 minutes. During this time you can put on the plaster **H**.

INTERPRETATION

NEGATIVE

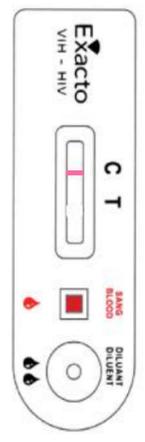POSITIVE

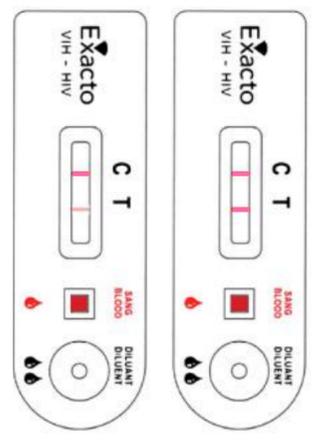PLACE YOUR TEST HERE

HIV-NEGATIVE LIKELY

If absence of Test band T and presence only of the Control band C (that indicates that the test was conducted properly).

- A **NEGATIVE** result can appear in a person who was infected over the last three months.
- You cannot conclude to be HIV-negative if you identify a risky exposure in the 3 months preceding this HIV screening self-test.

HIV-POSITIVE LIKELY

If presence of Test band T (more or less intense), and presence of the Control band C (that indicates that the test was conducted correctly).

- A **POSITIVE** result must be confirmed by a laboratory test. To do so, please contact your doctor or go to a screening centre.
- Pending confirmation of your result, use condoms to protect yourself and others.

- No band, or the test band T only : your EXACTO® TEST HIV is NOT VALID. Please do another test.
- Do not interpret the EXACTO Test HIV self-test if the blood sample has been transferred into the «DILUENT» round well and the diluent has been deposited in the «BLOOD» square well. In this case, the bottom of the «DILUENT» round well is red. You must repeat a new self-test or go to either a screening center or to your attending physician.
- Please refer to the national regulations or guidelines in force regarding the disposal of your test and its components. The disposable lancet may display an infectious risk after use. Therefore, the French legislation provides for its disposal in a needle box that is given free by your pharmacist. For more information, please ask your pharmacist for advice or visit www.exacto.fr or www.dastri.fr

PERFORMANCE

EXACTO® TEST HIV has undergone clinical assessments. A 99.9% correlation between **EXACTO® TEST HIV** and an ELISA laboratory test was observed. A clinical study* performed in hospitals using capillary blood showed that the **EXACTO® TEST HIV** has a sensitivity of 100% and a specificity of 100%. A practicability study conducted with a laymen population showed that 97.8% of participants obtained a correct and interpretable result. A readability study confirmed that 98% of the different types of results were interpreted correctly. **EXACTO® TEST HIV** can in rare clinical situations give erroneous results.

* HIV Self Test Exacto®, a Finger-Stick Whole Blood HIV1/2 Exhibit a very high Sensitivity and Specificity in both HIV infected known patients and HIV negative individuals, 2016. CHR Orléans T.PRAZUCK et al.

PRECAUTIONS FOR USE

- Do not freeze. Keep in a dry place between 2°C and 30°C.
- Remove the test cassette from the sealed pouch and use it within an hour.
- Run the test in a place where the temperature is between 15°C and 30°C.
- This test should be used only from a blood sample and should be performed immediately after blood sample collection at the fingertip.
- Single use
- The bottom of the square well BLOOD must be completely red before adding the diluent in the round well DILUENT.
- Keep the self-test and all components out of reach of children; ingesting the diluent may be dangerous and the lancet could potentially cause an injury.
- This test is designed for the self-screening of users who are not aware of their serological status and must not be used for the follow-up of antiretroviral treatments (ART) of HIV-positive persons. People using PrEP may obtain false negative results when using EXACTO® test HIV. A medical follow up with regular testing is preconized in this case.
- Negative results may be obtained in people infected less than three months before the test or infected with rare HIV variants. If you think you experienced a risky exposure over the last 3 months, you must visit your General Practitioner (GP) to get a prescription for a 4th generation ELISA-test in a laboratory, or use this self-screening test at least 3 months after the last risky exposure. Over this 3-month period, use condoms to protect yourself and others.
- A positive result must absolutely be confirmed by a laboratory test. Consult your doctor.
- The diagnosis of an infection by the HIV virus can only be done by a doctor on the basis of biological tests run in a clinical laboratory. **EXACTO® TEST HIV** is merely an extra tool aimed at encouraging the screening of HIV infections.
- **EXACTO® TEST HIV** cannot detect other infections (sexually transmitted infections or others). Erroneous results can be obtained if :
 - the instructions for use were not properly followed,
 - the aluminum pouch was damaged or the test was not run immediately after the pouch was opened,
 - storage conditions were not respected or the test was conducted after the expiry date indicated on the aluminum pouch.

Consult instructions for use	Sufficient for	Manufacturer
IVD in vitro diagnostic medical device	Use by	Do not reuse
Temperature limitation	LOT Batch code	REF Catalogue reference
DIL Diluent	Green point	TRIMAN unified signage

Date of last revision : 10/2019